

ONE-CARE™

SAFETY PEN NEEDLES



Insulin delivery has never been easier...

Safety

Clear safety shield automatically locks after use to protect the needle

Ergonomic Design

Ergonomic main housing with ribs for comfortable grip

Simple

Color indicator becomes invisible after engagement of safety mechanism

Universal

Compatible fit with most pen injectors

SAFE

- Automated shielding mechanism ensuring that needle is covered before and after use
- Single-use mechanism preventing sharps injuries and reuse of the device
- All safety pen needles are ETO-sterilized and protected with cover and sterility peel-off tab

SIMPLE

- Safety mechanism engagement indicator via "Click" sound and disappearance of color coded inner shield after fully inserting the needle into the skin
- Optimally designed, clear shield for precise placement on injection site with improved visibility

GENTLE

- Design with concealed needle helps to relieve patient's anxiety
- Polished, silicone-coated, ultra-thin needles are ideally positioned during skin penetration thanks to precise tracks for the safety mechanism elements enhancing patients comfort
- Ribbed housing provides a comfortable grip for easy attachment and removal from the pen

UNIVERSAL

- Universal design for compatible use with most injection pens available on the market
- Tight fit prevents leakages of insulin
- Easy to recognize color-coded product versions address a variety of patient skin types and injection requirements

UNIVERSAL STANDARDS

- Rigorously tested per international standards for fit, functionality and reliable performance

Compatibility


FlexTouch® Novo Nordisk
FlexPen® Novo Nordisk
NovoPen® Echo Novo Nordisk
NovoPen® 4 Novo Nordisk
NovoPen® 5 Novo Nordisk
Lantus SoloStar®
OptiPen® Pro 1 Sanofi
KwikPen® Humalog

Liprolog® Eli Lilly
Huminsulin® Eli Lilly
HumaPen® Luxura™ Eli Lilly
HumaPen® Luxura™ HD Eli Lilly
HumaPen® MEMOIR™ Eli Lilly
Exenatide Pen Byetta
Autopen® Classic Owen Mumford
Autopen® 24 Owen Mumford

Omnican® Pen 31B. Braun
Ypsopen® Ypsomed
OptiClik® Ypsomed
BerlinPen® Areo Berlin-Chemie
BerlinPen® Areo 2 Berlin Chemie
BerlinPen® 302 Berlin Chemie
Tactipen® Haselmeier
and many others



INSTRUCTIONS FOR USE

- Intuitive handling with simple steps. Insulin delivery has never been easier...

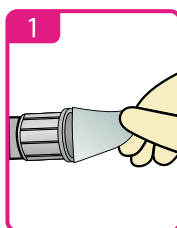



ATTACH

Peel off the sterility protector.
Attach the pen needle onto the pen by
screwing it in clockwise.
Remove the outer cover.

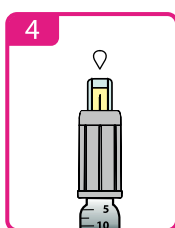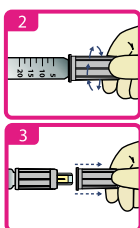

INJECT

Insert the needle into the skin at a 90-degree angle
in one continuous motion until clear shield retracts
and inject insulin by depressing the pen button.
The safety mechanism is triggered when the pen
needle is full pressed against the skin.

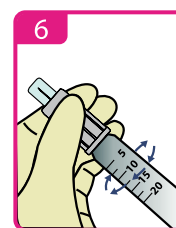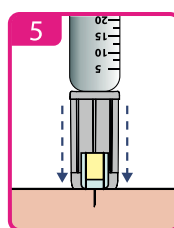

REMOVE AND DISPOSE

Lift the pen straight away from the skin.
Inner colored shield locks in place and becomes
invisible while the clear shield will automatically
cover the needle to protect users from injury.
Remove the needle from the pen by unscrewing it
counter-clockwise and simply dispose in a sharps
container.



ONE-CARE® Pen Needle with
sterility tab and outer cover



ONE-CARE® Pen Needle in the
ready for injection position



ONE-CARE® Pen Needle after use
(color-coded inner shield invisible)

PRODUCT SELECTION CHART

Available in two color-coded product versions for the right fit for every patient and skin type.

CODE	VERSION	DIAMETER (G) x DEPTH (mm/in)	PACKAGING
3101	Micro	30G x 5mm (3/16 in)	100 pcs / box 2,400 pcs / case
3102	Normal	30G x 8mm (5/16 in)	100 pcs / box 2,400 pcs / case

WWW.ONE-CARE.COM

