

ONE-CARE[®]

SAFETY PEN NEEDLES **PRO**



Insulin delivery has never been easier...

Safety

Clear safety shield automatically locks after use to protect the needle

Ergonomic Design

Ergonomic main housing with ribs for comfortable grip

Simple

Color indicator becomes invisible after engagement of safety mechanism

Universal

Compatible fit with most pen injectors

SAFE

- Automated shielding mechanism ensuring that needle is covered before and after use
- Single-use mechanism preventing sharps injuries and reuse of the device
- All safety pen needles are ETO-sterilized and protected with cover and sterility peel-off tab

SIMPLE

- Safety mechanism engagement indicator via "Click" sound and disappearance of color coded inner shield after fully inserting the needle into the skin
- Optimally designed, clear shield for precise placement on injection site with improved visibility

GENTLE

- Design with concealed needle helps to relieve patient's anxiety
- Polished, silicone-coated, ultra-thin needles are ideally positioned during skin penetration thanks to precise tracks for the safety mechanism elements enhancing patients comfort
- Ribbed housing provides a comfortable grip for easy attachment and removal from the pen

UNIVERSAL

- Universal design for compatible use with most injection pens available on the market
- Tight fit prevents leakages of insulin
- Easy to recognize color-coded product versions address a variety of patient skin types and injection requirements



INTERNATIONAL STANDARDS

- Rigorously tested per international standards for fit, functionality and reliable performance


Compatible with all type A insulin pens as per ISO11608-2, including:

COMPATIBILITY	Novo Nordisk	NovoPen® FlexPen®, NovoPen® FlexTouch®, NovoPen® Echo, NovoPen® 4, NovoPen® 5, InDuo®, InnoLet®, Victoza®
	Eli Lilly	KwikPen®, Forteo®, Huminsulin®, HumaPen®, HumaPen® Savvio™, HumaPen® Luxura™ HD, HumaPen® MEMOIR™
	Sanofi	Apidra® SoloSTAR®, ClickSTAR®, Lantus® SoloSTAR®, Toujeo® SoloSTAR®, AllSTAR™, JuniorSTAR®, Lyxumia®, OptiPen® Pro 1
	BD	BD® Pen
	Astra Zeneca	Byetta®
	Owen Mumford	Autopen® Classic, Autopen® 24
	B. Braun	Omnican® Pen 31
	Ypsomed	Ypsopen®, OptiClik®
	Berlin Chemie	BerliPen® Areo 2/3, BerliPen® Areo 301/302, BerliPen® Junior, BerliPen® Precision
	Haselmeier	Diapen®, TactiPen®
Bioton	GensuPen®	

Brands are trademarks of their respective owners

INSTRUCTIONS FOR USE

- Intuitive handling with simple steps. Insulin delivery has never been easier...

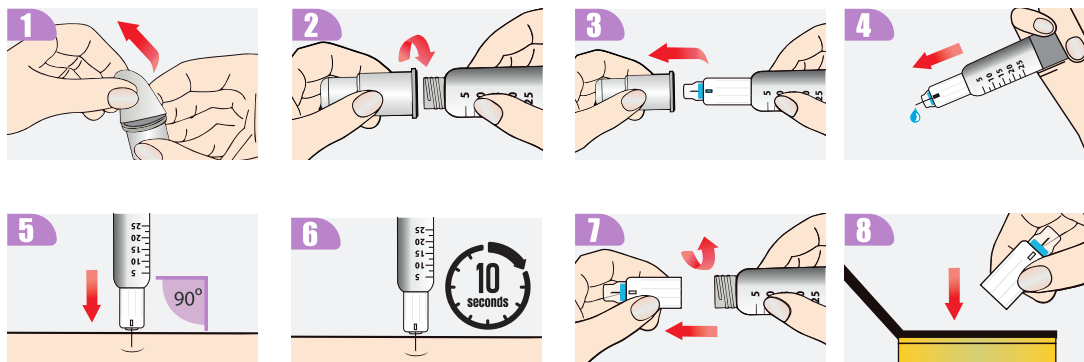

ATTACH

- 1 - Peel off the sterility seal
- 2 - Attach the pen needle onto the pen by screwing it in clockwise
- 3 - Remove the outer cover

PRIME AND INJECT

- 4 - Prime the pen needle by dispensing a drop of insulin
- 5 - Insert pen needle into injection site under 90 degree angle
- 6 - Wait 10 seconds after fully depressing pen button to deliver the entire insulin dose

REMOVE AND DISPOSE

- 7 - Remove the needle from the pen by unscrewing it counter-clockwise
- 8 - Dispose of the used pen needle into sharps container



PRODUCT SELECTION CHART

- Available in two color-coded product versions for the right fit for every patient and skin type

CODE	DESCRIPTION	DIAMETER (G) x DEPTH (mm/in)	COLOR	PACKAGING
3201	ONE-CARE PRO SAFETY PEN NEEDLES	31G x 5mm (3/16 in)	PURPLE	100 pcs / box 2,000 pcs / case
3202	ONE-CARE PRO SAFETY PEN NEEDLES	31G x 8mm (5/16 in)	BLUE	100 pcs / box 2,000 pcs / case

WWW.ONE-CARE.COM

